

Drugs of Abuse Detection		All products have a 25T/Box Pack Size		
Cat. No.	Description	Specimen	Format	CE Mark
D-001	Amphetamine	Urine	Dip Test/ Cassette	✓
D-001D	Amphetamine	Urine	Dip Test	✓
D-001S	Amphetamine	Saliva	Oral Collection Cassette	✓
D-002S	Methamphetamine	Urine	Dip Test/ Cassette	✓
D-002D	Methamphetamine	Urine	Dip Test	✓
D-002S	Methamphetamine	Saliva	Oral Collection Cassette	✓
D-003	Cocaine	Urine	Dip Test/ Cassette	✓
D-003D	Cocaine	Urine	Dip Test	✓
D-003S	Cocaine	Saliva	Oral Collection Cassette	✓
D-004	Morphine	Urine	Dip Test/ Cassette	✓
D-004D	Morphine	Urine	Dip Test	✓
D-004S	Morphine	Saliva	Oral Collection Cassette	✓
D-005	THC (Marijuana)	Urine	Dip Test/ Cassette	✓
D-005D	THC (Marijuana)	Urine	Dip Test	✓
D-005S	*THC (Marijuana) 20ng/ml	Saliva	Oral Collection Cassette	✓
D-006	Benzodiazepines	Urine	Dip Test/ Cassette	✓
D-006D	Benzodiazepines	Urine	Dip Test	✓
D-006S	Benzodiazepines	Saliva	Oral Collection Cassette	✓
D-007	PCP	Urine	Dip Test/ Cassette	✓
D-007D	PCP	Urine	Dip Test	✓
D-007S	PCP	Saliva	Oral Collection Cassette	✓
D-008	Barbiturates	Urine	Dip Test/ Cassette	✓
D-008D	Barbiturates	Urine	Dip Test	✓
D-008S	Barbiturates	Saliva	Oral Collection Cassette	✓
D-009	MDMA (Ecstasy)	Urine	Dip Test/ Cassette	✓
D-009D	MDMA (Ecstasy)	Urine	Dip Test	✓
D-009S	MDMA (Ecstasy)	Saliva	Oral Collection Cassette	✓
D-010	Oxycodone(OXY)	Urine	Dip Test/ Cassette	✓
D-010D	Oxycodone(OXY)	Urine	Dip Test	✓
D-010S	Oxycodone(OXY)	Saliva	Oral Collection Cassette	✓
D-011	TCA	Urine	Dip Test/ Cassette	✓
D-011D	TCA	Urine	Dip Test	✓
D-011S	TCA	Saliva	Oral Collection Cassette	✓
D-012	Methadone (MTD)	Urine	Dip Test/ Cassette	✓
D-012D	Methadone (MTD)	Urine	Dip Test	✓
D-012S	Methadone (MTD)	Saliva	Oral Collection Cassette	✓
D-013	Cotinine	Urine	Dip Test/ Cassette	✓
D-013D	Cotinine	Urine	Dip Test	✓
D-013S	Cotinine	Saliva	Oral Collection Cassette	✓
D-014	BUP	Urine	Dip Test/Cassette	
D-014S	BUP	Saliva	Oral Collection Cassette	
D-015	Syn-THC Synthetic	Urine	Dip Test/Cassette	

*THC TESTING AT PARENT Δ 9 at 20ng/ml

DRUGS OF ABUSE DETECTION

AMPHETAMINE (URINE/SALIVA)

This test can detect the presence of Amphetamines in the urine of a subject. Amphetamines generally stays in the system 13 days after use. The Cutoff concentration for this test is 1000 ng/ml, as recommended by the Substance Abuse and Mental Health Services Administration.

METHAMPHETAMINE (URINE)

This test is specifically designed to determine whether someone is using Methamphetamine vs. Amphetamines. The presence of methamphetamines can be detected up to 35 days after use.

COCAINE (URINE) / COCAINE (SALIVA)


The cocaine urine or saliva test is a standard single panel drug test which looks specifically for cocaine and its metabolite benzoylecgonine. The GC/MS cutoff level is 300 ng/ml and can be detected up to 23 days after use.

MORPHINE (URINE) / MORPHINE (SALIVA)

The Morphine urine or saliva test is a standard single panel drug test which looks specifically for opiate and opiate metabolites in either urine or saliva. Certain opiates can stay in the urine for up to 14 days.

THC 'MARIJUANA' (URINE) / THC (SALIVA)

This is a one step rapid, qualitative immunoassay for the detection of tetrahydrocannabinol compounds (9carboxyTHC) in the urine or saliva.

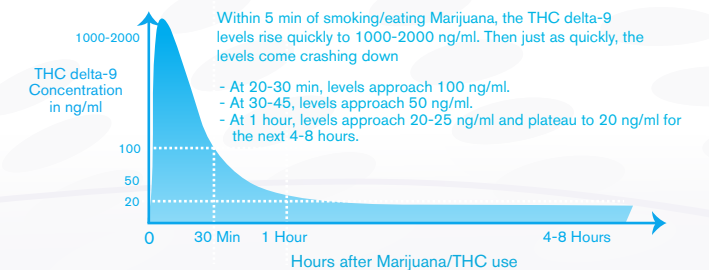

BENZODIAZEPINES (URINE) / BENZODIAZEPINES (SALIVA)

Benzodiazepines are a class of gabaminergic psychoactive chemicals used in the treatment of anxiety and sleep disorders. This class includes alprazolam (Xanax), diazepam (Valium), clonazepam (Klonopin) and others. Benzodiazepines are used extensively in emergency medicine and can be detected for 24 days after inception.

Continued Drugs of Abuse Detection		All products have a 25T/Box Pack Size		
Cat. No.	Description	Specimen	Format	CE Mark
D-015 S	Syn-THC Synthetic	Saliva	Oral Collection Cassette	
D-016	Mephedrone	Urine	Dip Test/Cassette	
D-016 S	Mephedrone	Saliva	Oral Collection Cassette	
D-017S	Alcohol	Saliva	Oral Collection Cassette	
D-020	2 in 1	Urine	Cassette	
D-021	3 in 1	Urine	Cassette	✓
D-022	4 in 1	Urine	Cassette	✓
D-023	5 in 1	Urine	Cassette	✓
D-024	6 in 1	Urine	Cassette	✓
D-025	7 in 1	Urine	Cassette	✓
D-026	8 in 1	Urine	Cassette	✓
D-027	9 in 1	Urine	Cassette	✓
D-028	10 in 1	Urine	Cassette	✓
D-029	11 in 1	Urine	Cassette	✓
D-030	12 in 1	Urine	Cassette	✓

*THC TESTING AT PARENT Δ 9 at 20ng/ml

PCP (URINE) / PCP (SALIVA)

PCP is a powerful psychedelic known for its dissociative effects at higher doses. It is best known for stories of the strange and sometimes violent behavior of people under its influence. PCP can stay in your system anywhere between 1 and 7 days. This test has a 25ng/ml cut off with a 17 day detection period.

BARBITURATES (URINE) / BARBITURATES (SALIVA)

Barbiturates are a class of sedative drugs prescribed for anxiety and insomnia. Barbiturates can stay in your system up to 3 days after use. This test can be completed in 5 min. with a cut off of 300 ng/l.

MDMA 'ECSTASY' (URINE) / MDMA 'ECSTASY' (SALIVA)

MDMA is one of the most popular recreational psychoactives, most commonly sold in the form of "ecstasy" tablets. It is known for its empathogenic, euphoric, and stimulant effects. It has also been used in psychotherapy. The substance can be detected 13 days after use. These tests look for the MDMA compound and can be sold as either a Urine or Saliva test.

OXYCODONE (URINE) / OXYCODONE (SALIVA)

Oxycodone is a semisynthetic opioid analgesic, used primarily in the treatment of pain. It saw a significant increase in recreational use when the higherdose timerelease form OxyContin was introduced in 1995. Some common types are: Percocet, Oxycontin, and Percodan. This test has a 100 ng/ml cut off in 5 minutes with a 34 day detection period.

TCA 'TRICYCLIC ANTIDEPRESSANTS' (URINE) / TCA (SALIVA)

Tricyclic antidepressants are the most common cause of death from prescription drugs. They can stay in your system for up to 10 days. Some common TCAs are: amitriptyline (Elavil), desipramine (Norpramin), imipramine (Tofranil), nortriptyline (Pamelor) and protriptyline (Vivactil). This test has a 1000 ng/ml cut off in 5 minutes with a ten day detection period.

METHADONE 'MTD' (URINE) / METHADONE 'MTD' (SALIVA)

A Methadone specific test must be performed to detect Methadone. Methadone is used as a substitute for opiates. The detection period in urine is 1-10 days with a cutoff level of 300 ng/ml.

COTININE (URINE) / COTININE (SALIVA)

Cotinine is a biomarker for exposure to tobacco smoke and is sometimes used as an antidepressant. This test has a 200 ng/ml cutoff in 5 minutes with a 7 day detection period.

BUPRENORPHINE 'BUP' (URINE) / BUPRENORPHINE 'BUP' (SALIVA)

Buprenorphine is a semisynthetic opioid derived from thebaine. It is an analgesic (pain killer) and has gained a great deal of notoriety for its ability to interrupt severe opiate addiction, including heroin and methadone addictions. This test has up to a 10 dya detection period.

SYN-THC SYNTHETIC (URINE) / SYN-THC SYNTHETIC (SALIVA)

Synthetic cannabis is a psychoactive designer drug derived from natural herbs sprayed with synthetic chemicals that, when consumed, allegedly mimic the effects of cannabis. This drug stays in your system up to 3 days after use, and has commonly been used by those attempting to bypass standard drug testing.

MEPHEDRONE (URINE) / MEPHEDRONE (SALIVA)

4-Methylmethcathinone is a synthetic stimulant with empathogenic effects. Between 2007 and 2009, it became available for purchase online, was used increasingly in Europe, Australia, and New Zealand, associated with several deaths, and, as a result, was controlled in some countries. Because of the similarity in names, it is sometimes confused with 'Methedrone' (4-methoxymethcathinone) or 'Methylone' (bkMDMA). Mephedrone is a stimulant substitute similar to amphetamines or MDMA, it can stay in your system for 27 days.

"27 day detection period, stimulant substitute."

ALCOHOL (SALIVA)

Alcohol has a 6-24 hours detection period in urine. This test can also detect Ethyl glucuronide (EtG), a direct metabolite of alcohol, which allows an extended window for detecting consumption.

2 IN 1 CASSETTE, 3 IN 1 CASSETTE, 4 IN 1 CASSETTE, 5 IN 1 CASSETTE,

6 IN 1 CASSETTE, 7 IN 1 CASSETTE, 8 IN 1 CASSETTE, 9 IN 1 CASSETTE,

10 IN 1 CASSETTE, 11 IN 1 CASSETTE, 12 IN 1 CASSETTE,



DRUG DETECTION CUPS

Drug Detection Cups		All products have FDA (510k) and CE Mark		
Cat. No.	Description	Specimen	Format	Drugs Tested
D-104	4 Panel Drug Cup	Urine	Drug Cup	4
D-105	5 Panel Drug Cup	Urine	Drug Cup	5
D-106	6 Panel Drug Cup	Urine	Drug Cup	6
D-107	7 Panel Drug Cup	Urine	Drug Cup	7
D-108	8 Panel Drug Cup	Urine	Drug Cup	8
D-109	9 Panel Drug Cup	Urine	Drug Cup	9
D-110	10 Panel Drug Cup	Urine	Drug Cup	10
D-111	11 Panel Drug Cup	Urine	Drug Cup	11
D-112	12 Panel Drug Cup	Urine	Drug Cup	12
D-112	12 Panel Drug Cup (CLIA Waived)	Urine	Drug Cup	12

A cup for urine testing containing any combination up to 12 panel drug cup tests of your choosing. Can come with a variety of Adulterants such as: Oxidants, Specific Gravity, pH, Nitrite, Glutaraldehyde, Creatinine

"Custom Combination, Adulterants: Oxidants, Specific Gravity, pH, Nitrite, Glutaraldehyde, Creatinine."

4 PANEL DRUG CUP (FDA 510K)

5 PANEL DRUG CUP (FDA 510K)

6 PANEL DRUG CUP (FDA 510K)

7 PANEL DRUG CUP (FDA 510K)

8 PANEL DRUG CUP (FDA 510K)

9 PANEL DRUG CUP (FDA 510K)

10 PANEL DRUG CUP (FDA 510K)


11 PANEL DRUG CUP (FDA 510K)

12 PANEL DRUG CUP (FDA 510K)


12 PANEL DRUG CUP (FDA 510K CLIA WAIVED)

QIKTECH INSTANT MULTI-DRUG DETECTION CUP

TEST INSTRUCTIONS

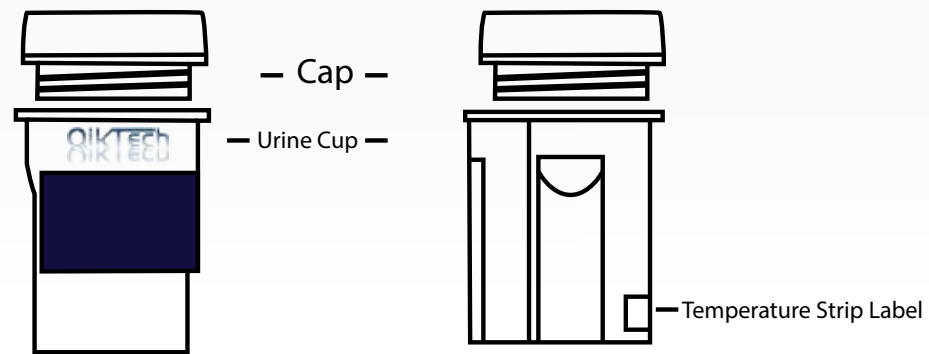


- 1 Remove Cap.
- 2 Donor provides specimen
- 3 Replace the Cap securely. Peel off the blue label that covers the result window.

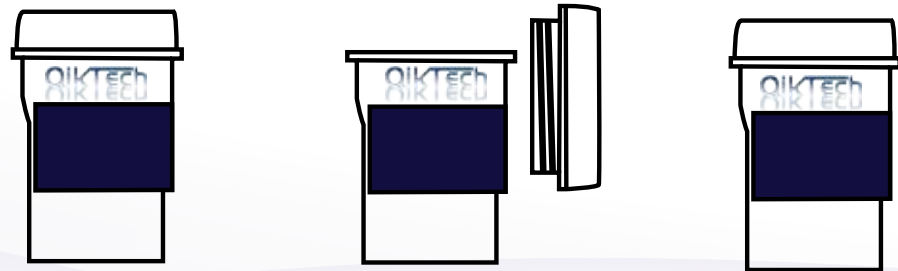

- 4 Read results at 5 minutes. Test results remain stable for a minimum of 60 minutes.



Negative results can be read as soon as they appear.
Read positive results at 5 minutes.

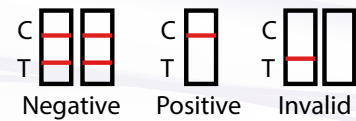

All positive results are presumptive and should be confirmed by an alternative method (e.g. GC/MS). A photocopy of the results can be made by placing the entire sealed cup with the card face down on the copier.